

# Estimation of $\text{Zn}^{2+}$ & $\text{Pb}^{2+}$ ions

---



IN A GIVEN SOLUTION COMPLEXOMETRICALLY

**ROHAN SINGH, NUTAN SHARMA**

Department of Chemistry

University of Delhi


# Chemicals Required

- ❖ Complexone: EDTA solution
- ❖ Standard: ZnO as primary standard
- ❖ Indicators: Eriochrome Black-T (EBT) and Xylenol Orange (XO)
- ❖ Buffers:  $\text{NH}_4\text{OH}$ - $\text{NH}_4\text{Cl}$  solution and Hexamine
- ❖ Auxiliaries: Ammonium acetate and acetic acid


# Complexometric Titrations

---

- ❖ Complexometric titrations are a form of volumetric analysis in which the formation of a colored complex is used to indicate the end point of a titration
- ❖ These are used mainly to determine metal ions in a mixture by use of complex-forming reactions
- ❖ Most common chelating agent used for in these titrations is EDTA

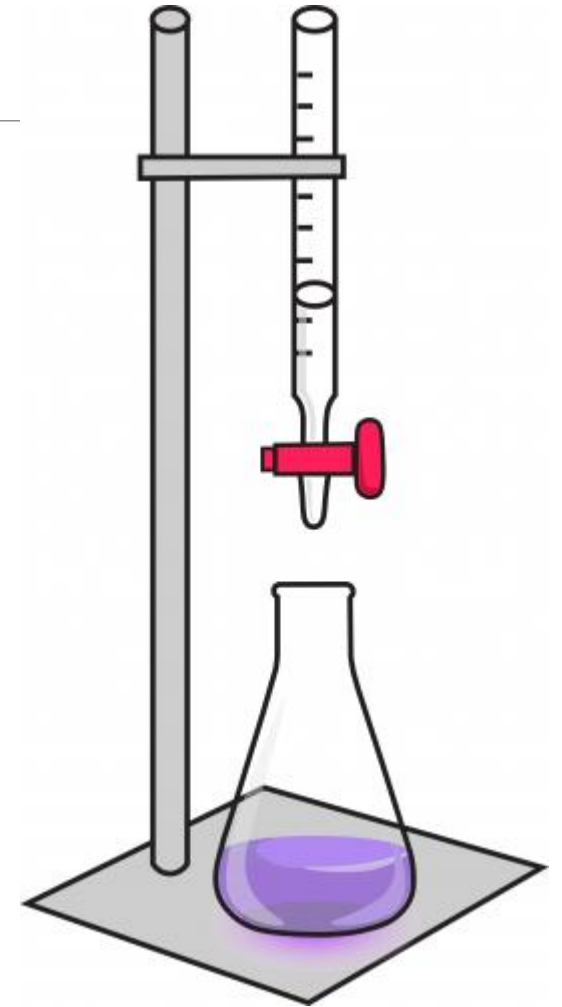

# Why EDTA??

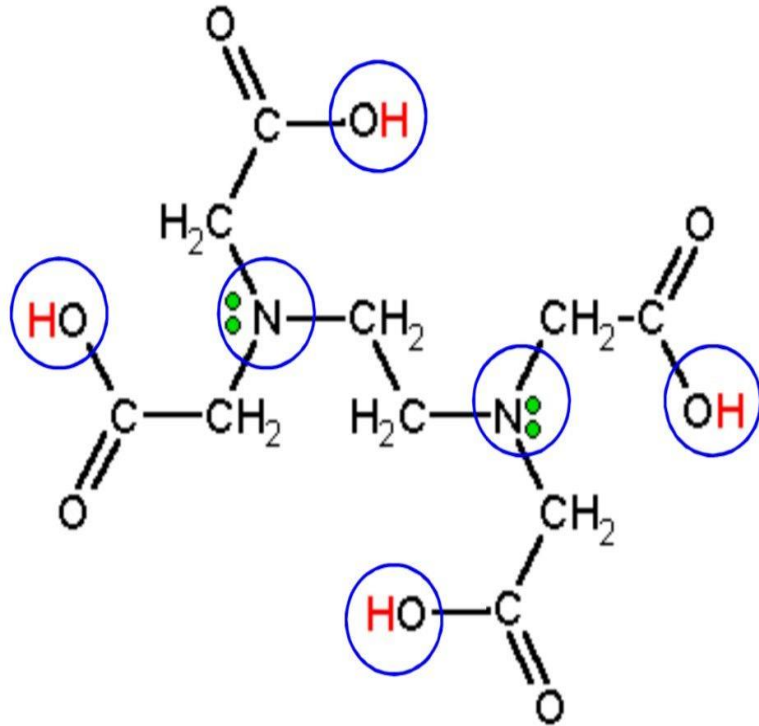

6 metal binding sites making it  
hexadentate ligand

- Reasons for using EDTA
  - Forms strong 1 : 1 complexes with most metal ions
  - Large formation constants with many metal ions
  - Forms stable, water-soluble metal complexes
- Titrations with EDTA gets affected by pH factor, stability constant of the complex formed etc.
- pH is important to guarantee the selectivity and accuracy of the analysis
- The more stable the chelate, the lower the pH at which the titration can be performed.

# Estimation of $\text{Zn}^{2+}$ & $\text{Pb}^{2+}$

IN A MIXTURE

- ❖  $\text{Zn}^{2+}$  and  $\text{Pb}^{2+}$  when present **simultaneously** in a solution can be estimated using **direct titration** with EBT and Xylenol Orange indicator in different pH ranges
- ❖ **Different pH range** is preferred for increasing selectivity
- ❖ Before proceeding to titration,  $\text{Zn}^{2+}$  and  $\text{Pb}^{2+}$  ions must be **separated** using a suitable precipitating agent.
- ❖ Another technique that might also be used is **masking** which eliminates the need for separating both the ions
- ❖ Cyanide ion can be used as masking agent in this method.

# Need for separation of $\text{Zn}^{2+}$ and $\text{Pb}^{2+}$

Formation constant of metal-EDTA complex:

$$K_f = \frac{[\text{M-EDTA}]}{[\text{M}][\text{EDTA}]}$$

A metal, M can be very easily titrated in presence of other metal, N if the formation constant values of these ions is such that:

$$K_M/K_N \geq 10^8 \quad (\text{in presence of complex forming indicator})$$

But

$$\begin{aligned} K_f(\text{Zn-EDTA complex}) &= 10^{16.5} \\ K_f(\text{Pb-EDTA complex}) &= 10^{18.0} \end{aligned}$$

$$\text{making } K_f(\text{Pb})/K_f(\text{Zn}) = 10^{1.5}$$

Because of this **proximity of formation constant values** of these two complexes,  $\text{Zn}^{2+}$  and  $\text{Pb}^{2+}$  must be separated before proceeding to titration. (formation of both Zn-EDTA and Pb-EDTA complexes in the same solution makes it difficult to estimate their concentrations individually)

# Separation of $\text{Pb}^{+2}$

- Precipitating agent is chosen in such a way that it forms a precipitate with only one of the metal ion leaving out the other in the solution form
  - **Precipitating agent:  $\text{K}_2\text{SO}_4$  solution.** In case  $\text{K}_2\text{SO}_4$  is unavailable,  $\text{K}_2\text{SO}_3$  or  $\text{Na}_2\text{SO}_3$  can be used after oxidation with  $\text{HNO}_3$
  - Conditions for precipitation:
    - ionic product > solubility product ( $K_{\text{sp}}$ )
    - Size Compatibility between cation and anion
    - lattice energy > hydration energy
- $\text{PbSO}_4$  satisfies all the 3 criteria**

$K_{\text{sp}} = [\text{Pb}^{2+}] \times [\text{SO}_4^{2-}] = 2.53 \times 10^{-8} \text{ M}^2$

Size of  $\text{Pb}^{2+} = 119 \text{ pm}$

Solubility = 0.00443 g/100 mL (20 °C)
- $\text{Pb}^{2+} + \text{SO}_4^{2-} \rightarrow \text{PbSO}_4$  : much more stable, white ppt is obtained
  - Water's dipole strength is too weak to pull away the ions (both anions and cations) from the strong crystals of lead sulphate
  - Solution is warmed to increase rate of precipitation.  $\text{K}_2\text{SO}_4$  is added dropwise to avoid super saturation
  - $\text{ZnSO}_4$  if formed, is soluble in water

# Digestion of precipitate


Occurs when a freshly-formed precipitate is left, usually at a higher temperature, in the solution from which it is precipitated.

**Need:** Results in bigger particles, generally purer, easier to wash & filter → **Ostwald ripening**

larger particles → greater volume to surface area ratio → lower energy state (more favored)

Molecules on the surface of a particle → energetically less stable than the ones already well ordered and packed in the interior.

System tries to lower its overall energy → molecules on the surface of a small particle tend to diffuse through solution and add to the surface of larger particle.

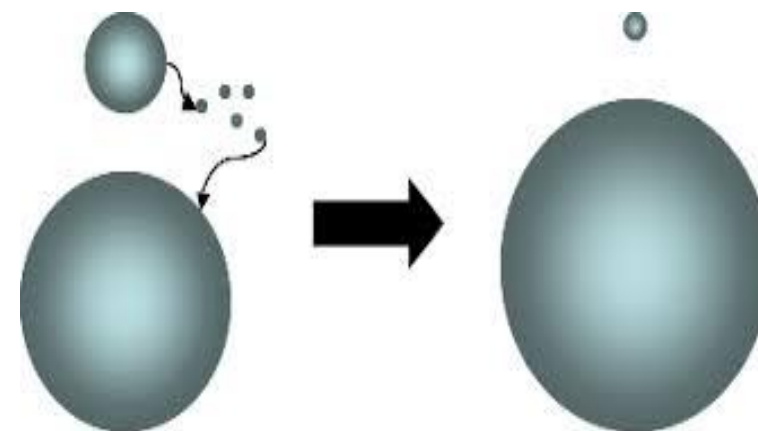



Energetic factors → large precipitates grow, drawing material from the smaller precipitates, which shrink



# Estimation of $\text{Zn}^{2+}$ ions

- $\text{Zn}^{2+}$  ions in the solution are estimated by titrating the filtrate with EDTA solution using EBT as indicator.
- $K_f (\text{Zn-EBT}) < K_f (\text{Zn-EDTA})$  { Difference in number of ligating centres in EBT and EDTA
- pH range = 10 ( $\text{NH}_4\text{OH-NH}_4\text{Cl}$  buffer is added)
  - For neutralizing effect of  $\text{H}^+$  ions released on reaction of  $\text{Zn}^{2+}$  with EDTA and EBT: Enables the formation of the  $\text{Zn}^{2+}$  with EDTA at the appropriate pH
  - avoid the  $\text{Zn}^{2+}$  ion precipitation (as  $\text{Zn}(\text{OH})_2$ ) in alkaline pH

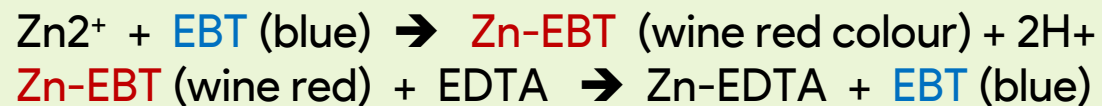


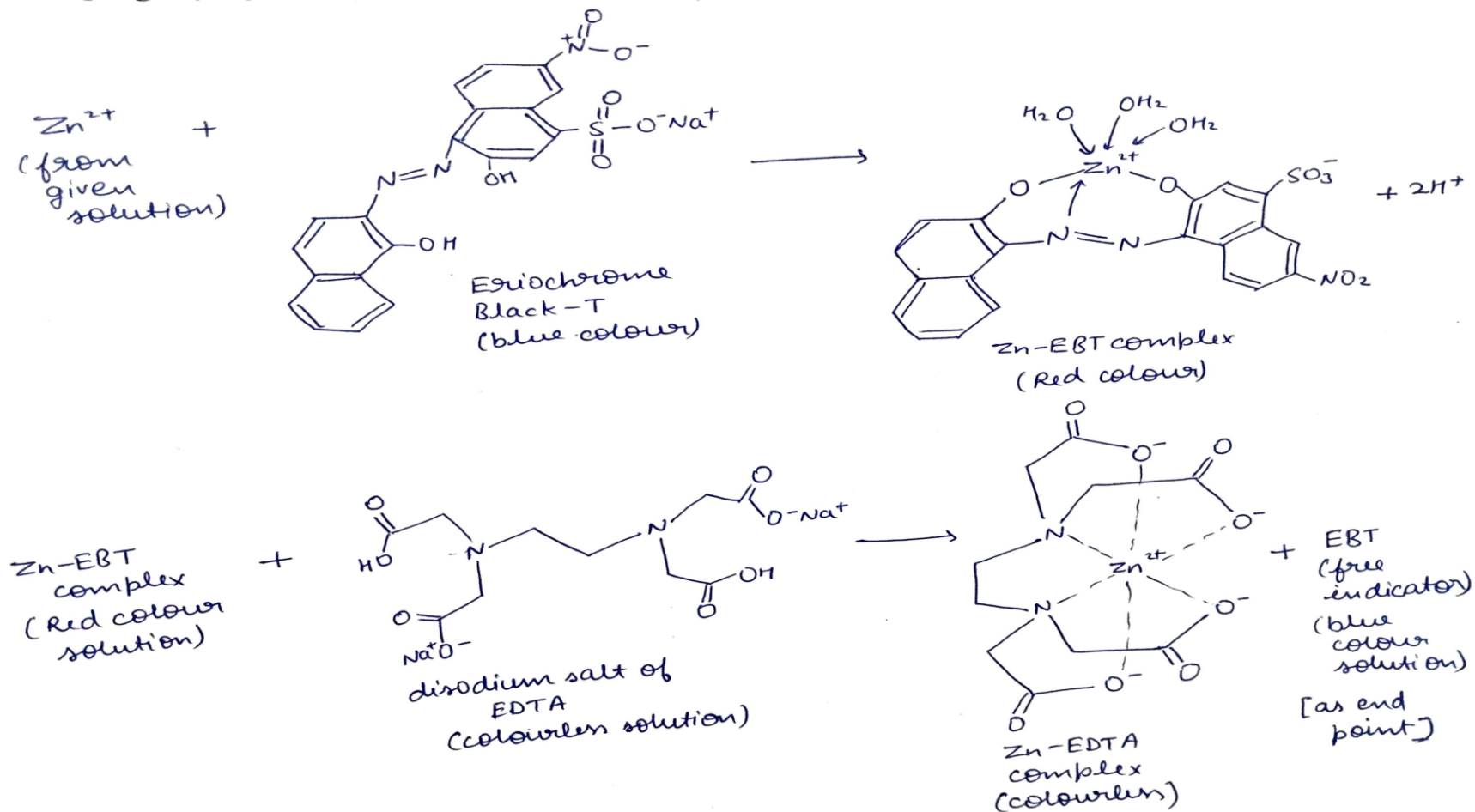


## 9. Reactions for estimation of $\text{Zn}^{2+}$


# Dissolution of $\text{PbSO}_4$ ppt

- Dissolved using ammonium acetate and acetic acid

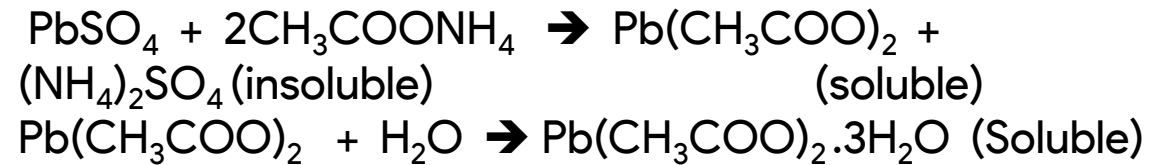


- Why lead acetate is soluble???
- It is soluble because of hydrogen bonding present due to its structure. In the trihydrate, the  $\text{Pb}^{2+}$  ion's coordination sphere consists of nine oxygen atoms belonging to 3 water molecules, two bidentate acetate groups and two bridging acetate groups.
- **Uses:** keeps  $\text{Pb}^{2+}$  in solution form preventing its hydrolysis

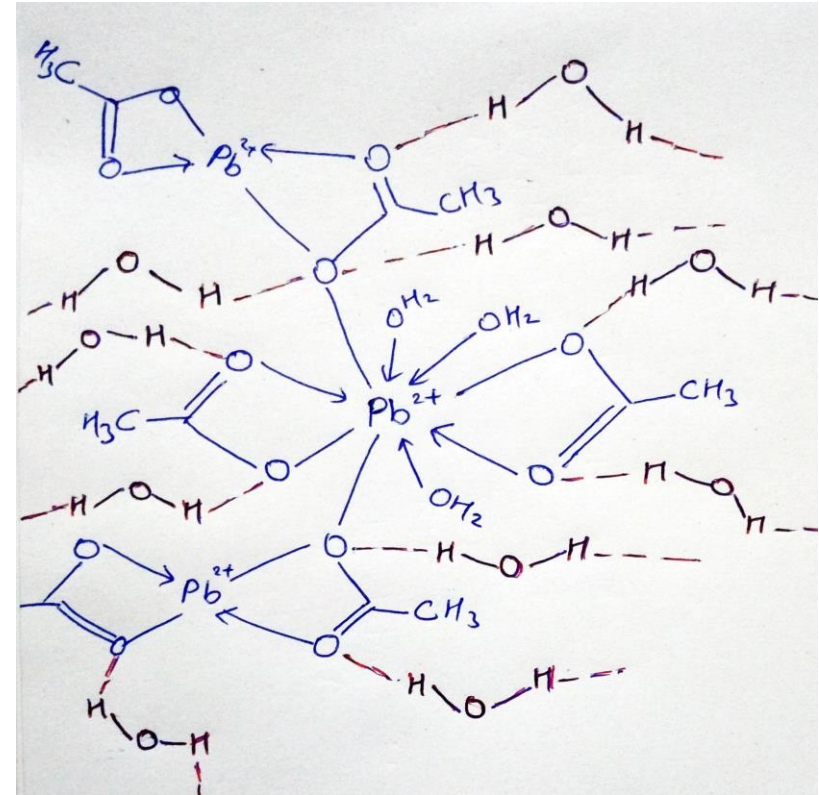

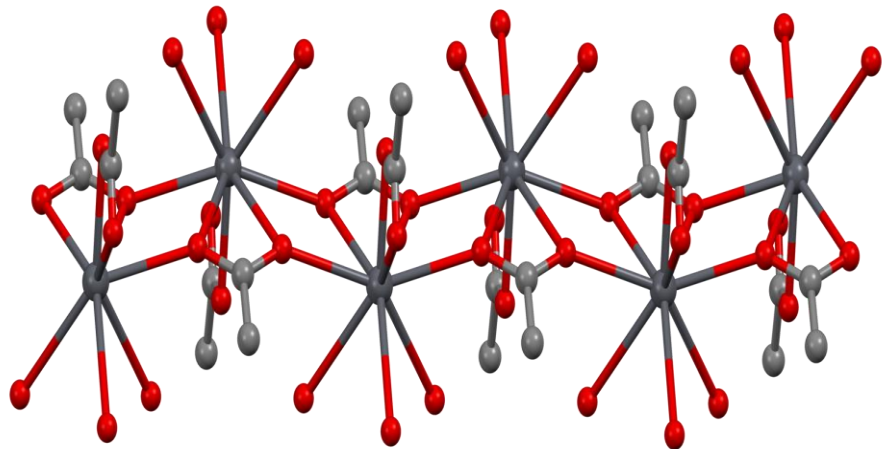


Hydrogen bonding in the crystal structure

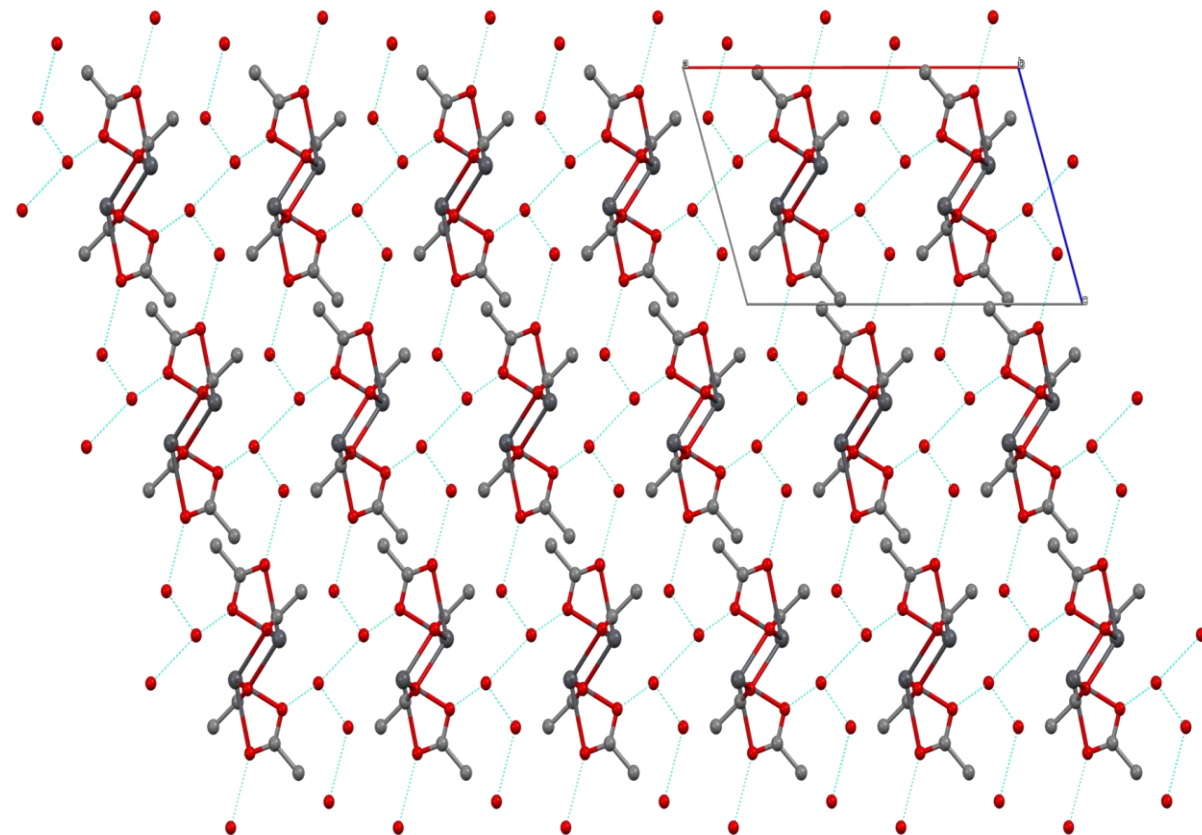


## 11. Structure of Lead Acetate (trihydrate)

# Estimation of $\text{Pb}^{2+}$ ions

- $\text{Pb}^{2+}$  ions in the solution can then be estimated by direct titration with EDTA solution using **Xylenol Orange** as indicator.
- pH maintained is 5  $\rightarrow$  obtained by adding hexamine buffer solution
- Why 5-6 pH??
  - For maximum stability of Pb-EDTA complex
  - to increase selectivity of EDTA
  - To prevent hydrolysis of  $\text{Pb}^{2+}$  ions

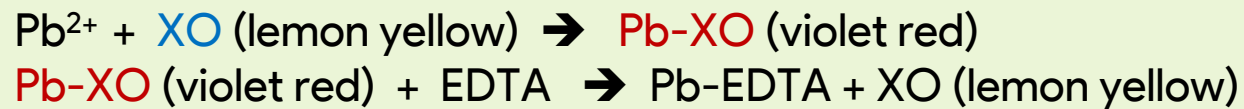

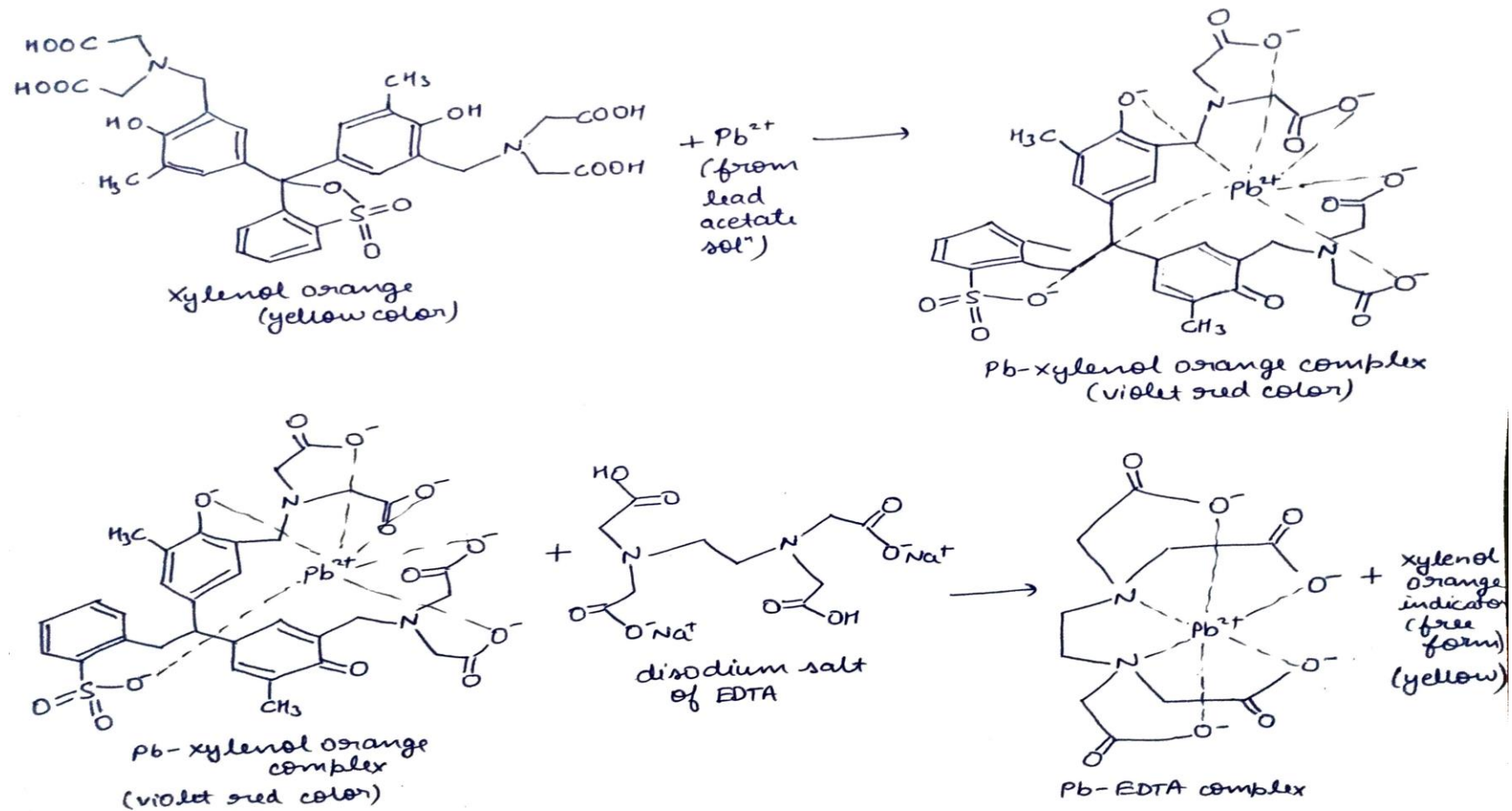

# 13. Reactions for estimation of $Pb^{2+}$




# Procedure

## 1. Preparation of standard solution of ZnO in 100 mL

Calculated amount of ZnO + Few drops of HCl (to make ZnO completely dissolve in water) + Distilled water up to 100 mL mark → Mixture is shaken well to get homogenous solution




## 2. Standardization of EDTA solution with std ZnO sol: EDTA solution is taken in burette.

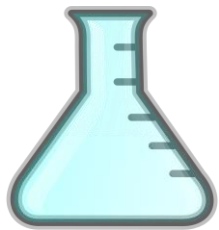

10 mL std ZnO sol + 5 mL of  $\text{NH}_4\text{OH-NH}_4\text{Cl}$  buffer



On adding a few drops EBT indicator, Zn-EBT complex (wine-red) is formed

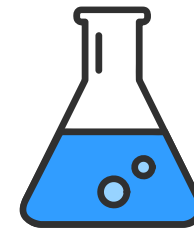

At end point, EDTA replaces EBT from Zn-EBT complex giving blue colour

### 3. Separation of $\text{Pb}^{2+}$ and $\text{Zn}^{2+}$ solution

- 20 mL of distilled water is added and heated up to  $60\text{ }^{\circ}\text{C}$ .
- 10%  $\text{K}_2\text{SO}_4$  solution is added dropwise until completion of precipitation
- Solution is digested for about 10 min
- A drop of  $\text{K}_2\text{SO}_4$  was added to check for completion of precipitation
- Solution is filtered using whatmann filter paper




Precipitate of  $\text{PbSO}_4$

$\text{Zn}^{2+}$  solution

### 4. Estimation of $\text{Zn}^{2+}$ solution


- Solution in burette: standardized EDTA

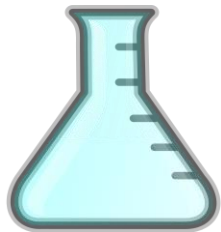

10 mL filtrate + 5 mL of  $\text{NH}_4\text{OH}$ - $\text{NH}_4\text{Cl}$  buffer



On adding a few drops EBT indicator, Zn-EBT complex (wine-red) is formed

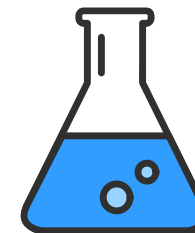

At end point, EDTA replaces EBT from Zn-EBT complex giving blue colour




5. Estimation of  $\text{Pb}^{2+}$  in given solution

- Precipitate of  $\text{PbSO}_4$  is dissolved using  $\text{CH}_3\text{COONH}_4$  &  $\text{CH}_3\text{COOH}$
- Solution is made up to mark of 100 mL
- Solution in burette: standardized EDTA solution



10 mL of  $\text{Pb}^{2+}$  sol. +  
hexamine buffer

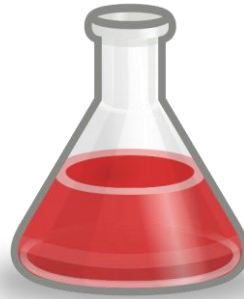

On adding a few drops of XO  
indicator,  $\text{Pb-XO}$  complex (violet-  
red) is formed



EDTA replaces XO from  $\text{Pb-XO}$   
complex  $\rightarrow$  free XO (lemon-  
yellow) is obtained

# Calculations

Let normality of std ZnO solution prepared be  $N_{\text{ZnO}}$ .

❖ **Normality of given EDTA solution:** *Zn-EBT (wine red) + EDTA → Zn-EDTA + EBT (blue)*

- $N_{\text{EDTA}} \times V_{\text{EDTA}} = N_{\text{ZnO}} \times V_{\text{ZnO}}$
- $N_{\text{EDTA}} = (N_{\text{ZnO}} \times 10) / V_{\text{EDTA}}$  ( $V_{\text{EDTA}}$  obtained from step 2)

❖ **Normality of  $\text{Zn}^{2+}$  in given solution:** *Zn-EBT (wine red) + EDTA → Zn-EDTA + EBT (blue)*

- $N_{\text{EDTA}} \times V_{\text{EDTA}} = N_{\text{Zn}^{2+}} \times V_{\text{Zn}^{2+}}$
- $N_{\text{Zn}^{2+}} = (N_{\text{EDTA}} \times V_{\text{EDTA}}) / 10$  ( $V_{\text{EDTA}}$  obtained from step 3)

❖ **Normality of  $\text{Pb}^{2+}$  in given solution:** *Pb-XO (violet red) + EDTA → Pb-EDTA + XO (yellow)*

- $N_{\text{EDTA}} \times V_{\text{EDTA}} = N_{\text{Pb}^{2+}} \times V_{\text{Pb}^{2+}}$
- $N_{\text{Pb}^{2+}} = (N_{\text{EDTA}} \times V_{\text{EDTA}}) / 10$  ( $V_{\text{EDTA}}$  obtained from step 4)


# Applications

OF TODAY'S EXPERIMENT

- ❖ Analysis of various lead-containing alloys such as bronzes, ounce metal, solders and tin- and lead-based alloys
- ❖ Analysis of content in zinc-lead ore
- ❖ Analysis of glass certified reference material– optical glasses etc.
- ❖ Removal of zinc and lead contamination from water, sewage and soil by estimating them

# Precautions

- ❖ Lead is a highly **poisonous metal** (whether inhaled or swallowed), affecting almost every organ and system in the human body. Hence masks and gloves are must while handling solution containing  $\text{Pb}^{2+}$  ions.
- ❖ **Cu, Co, Ni, Cr, Fe, Al** even in traces must be absent as they interfere in the titration by forming more stable M-EDTA complexes.
- ❖ Deionized water is preferred over distilled water.
- ❖  $\text{PbSO}_4$  ppt must be washed properly so as to make sure removal of  $\text{Zn}^{2+}$  solution from it.
- ❖ Burette must be rinsed with EDTA solution before titrations.

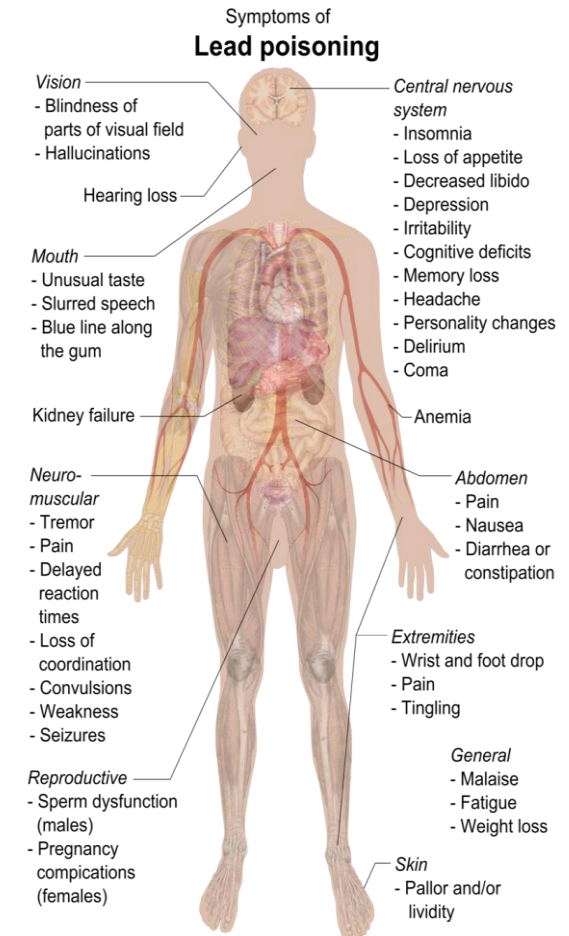

# References

- ❖ Vogel's textbook of quantitative chemical analysis
- ❖ Quantitative chemical analysis by Daniel C. Harris
- ❖ [https://en.wikipedia.org/wiki/Lead\(II\)\\_acetate](https://en.wikipedia.org/wiki/Lead(II)_acetate)
- ❖ <https://chemistry.stackexchange.com/>
- ❖ <https://www.chemeurope.com/en/encyclopedia/>